

Your Trustworthy ECG Assistants **SE-1200 Express**

12-Channel Electrocardiograph



SE-1200 Express

12-Channel Electrocardiograph



SE-1200 Express

SE-1200 Express (Basic)



Trolley



Rolling Stand



★ We still offer a basic version of SE-1200 Express to help with your budget limitation.



ECG Reports

Streamlined Workflow

Patient Information



User-friendly Design

- 8" high resolution color touch screen
- Alphanumeric keyboard

Streamlined Workflow

- One-touch operation with age and gender shortcuts
- Support linear/2D barcode scanner to simplify workflow

Flexible Printing Solutions

- Built-in high resolution thermal recorder to print on normal thermal paper and A4 size fax paper
- Support external USB printer

Data Management

- Internal storage up to 800 ECGs and can be enlarged by USB flash disk
- SE-1515 Data Management System software for data management on PC (optional)

Accurate Clinical Performances

- Weak signal detection with the advanced 24-bit A/D converter technology
- Anti-noise technology with comprehensive filters and high CMRR over 140dB
- ECG Signal Quality Indication
- 120s real-time ECG waveform freezing and review to assist diagnosis
- Automatic measurements and interpretation tested with CSE/AHA/MIT databases
- Extend printing once arrhythmia is detected for further rhythm analysis
- Built-in Glasgow Interpretation
- Advanced ECG analysis: Vector ECG (optional)

Extensive Communications

- DATA transmitting to PC via LAN/ internal Wi-Fi (Wi-Fi is optional)
- Report formats: PDF and optional SCP/FDA-XML/DICOM data export
- Bi-directional data communication with DMS

Seamless Connectivity


Network your ECG data with the help of EDAN's SE-1515 Data Management System in the same LAN. It is convenient to share, review, edit and diagnose ECG reports between different doctors. Meanwhile, with the help of Smart ECG Viewer, EDAN SE-1200 Express series are able to access EMR/PACs/GDT-based Systems via the standard protocols: DICOM Worklist/HL7/GDT.


Hospital System Access

Patient information (worklist) will be transmitted from hospital systems to ECG machines. Meanwhile, reports or final diagnostic results will be uploaded to hospital systems for archiving.

Standard Communication Protocols

HL7/DICOM worklist are supported on SE-1515 Data Management System.

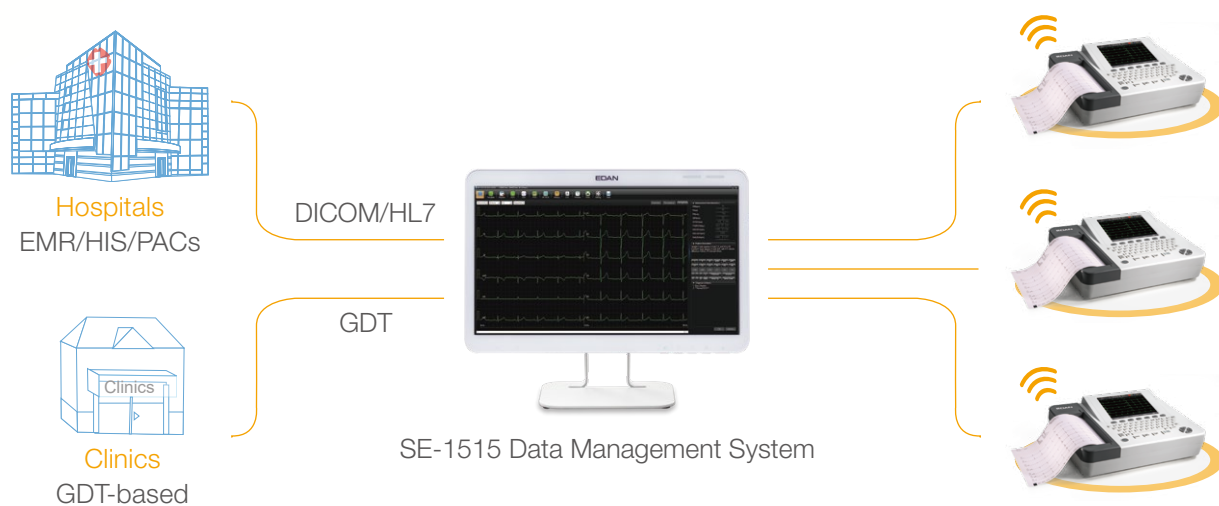

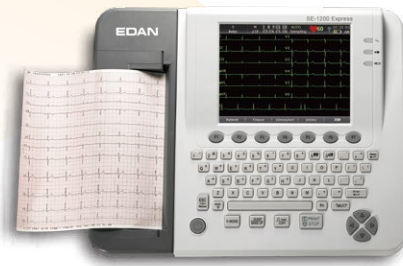

SE-1200 Express 12-Channel Electrocardiograph



Global Headquarters:

Edan Instruments, Inc. | 15 Jinhui Road, Pingshan District, Shenzhen
518122 P.R. China | +86.755.26898326 | www.edan.com | info@edan.com

U.S. and Canada inquiries:

EDAN Diagnostics, Inc. | 9918 Via Pasar, San Diego, CA 92126
+1.858.750.3066 | www.edandiagnostics.com | edan-info@edandiagnostics.com

© Edan Instruments, Inc. All rights reserved. Features and specifications are subject to change without prior notice. No reproduction, copy or transmission may be made without written permission. Not all products or features are available in all countries, contact Edan for local availability.



ENG-ECG-SE-1200 Express
V3.6-20240605