

Specialist of *Sleep* Diagnosis and Therapy



CE 0123

3.5-inch Color LCD Display



With an exquisite control knob, BMC RESmart Generation II Series bring you unsurpassed operational experience. Real time flow wave connects us to better patient care and compliance.

Eco Smart Technology

Eco Smart heating system with innovative dual water chambers design, Storage Chamber VS Heating Chamber. Optimum Humidity and Temperature are delivered to patient after traveling through the heating chamber. Accurate water quantity control, with real time compensation from Storage Chamber. Eco Smart warm up, high efficiency, lower water bills.


iCode II

Compliance Data is encoded by QR code
Google play and App Store apps are available
High reliability – recognition rate up to almost 100%
High performance – high recognition speed.
More information – less information input

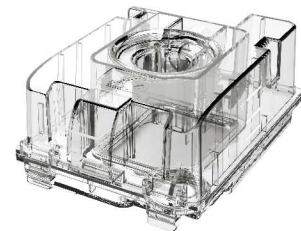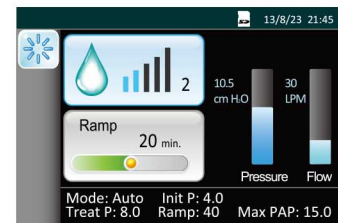

Modes	Parameter	RESmart GII CPAP	RESmart GII-Auto CPAP
	Cover Color	Black	Black
	LCD Screen	2.4 inches	3.5 inches
	Mode	CPAP	CPAP/Auto
	Pressure Range	4-20 cm H ₂ O	4-20 cm H ₂ O
	Pressure Increment	0.5 cm H ₂ O	0.5 cm H ₂ O
	Reslex	Off/1/2/3	Off/1/2/3
	Ramp	0-60 minutes	0-60 minutes
	Sound	<30 dBA	<30 dBA
	Alarms	Power Fail/ High Leak	Power Fail/ High Leak
	Humidity Level	Off/1-5	Off/1-5
	Auto On/Off	Off (Default)	Off (Default)
	Data Storage	SD Card	SD Card
	Dimensions (L×W×H)	6.6" x 7.7" x 4.6" (170 × 196 × 118 mm)	
	Weight	3.3 lb (1.5 kg) offset	
	Power Supply	Input range: 100-240V; 50-60Hz; Power consumption: <60 VA	

ECO Smart

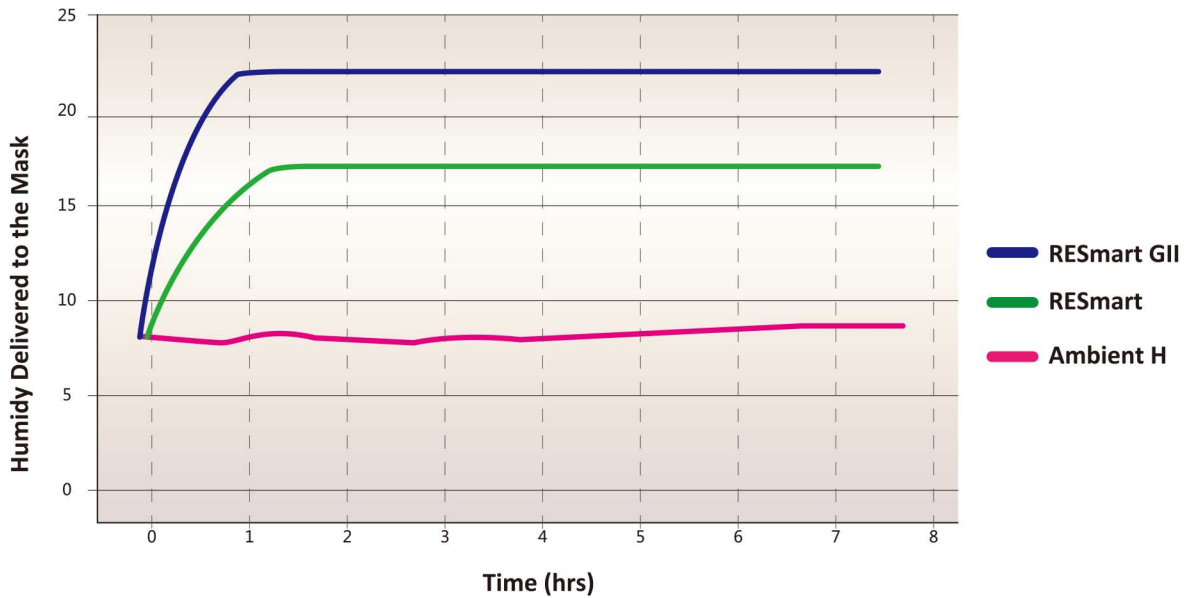

BMC Medical Co., Ltd.

Add: 5/F Main Building, No.19 Gucheng Street West,
Shijingshan, Beijing 100043, China
Tel: +86-10-51663880
Fax: +86-10-51663880 ext. 810
www.bmc-medical.com