

Building on more than 30 years of experience in aerosol therapy, Philips Home Nebulizer offers high quality aerosol therapy for the effective management of respiratory conditions, at an affordable price.

Reliable and easy to use, Home Nebulizer is suitable for nebulizing most aerosol therapy drugs.















Benefits



Economical choice

Affordable aerosol therapy – suitable for children and adults...



Effective therapy

Reliable nebulizer therapy delivered using the clinically proven SideStream nebulizer technology**.



Treatment in 6 mins*

A compact design that provides delivery of consistent aerosol medication treatments in as little as 6 minutes.



Easy to use

The smallest of the Philips home compressor nebulizers, Home Nebulizer is also lightweight and simple to operate.














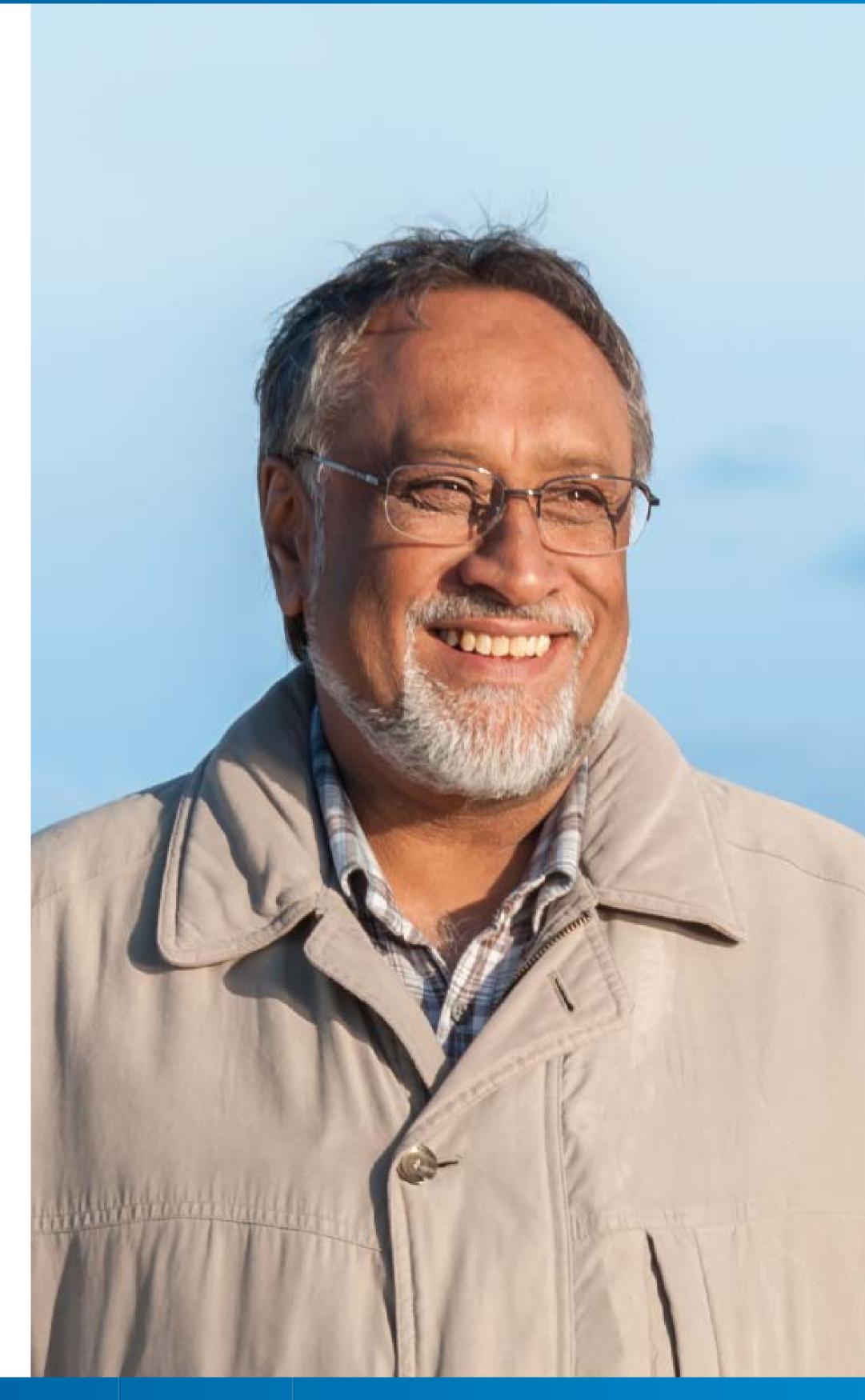
















Philips Home Nebulizer is combined with our highly efficient SideStream nebulizer. SideStream has been designed to draw more air in, via an active venturi system, resulting in fast drug delivery and short treatment times. It also has a unique diamond jet design which helps to minimize wear from frequent use and aids consistent delivery of a high-quality aerosol.

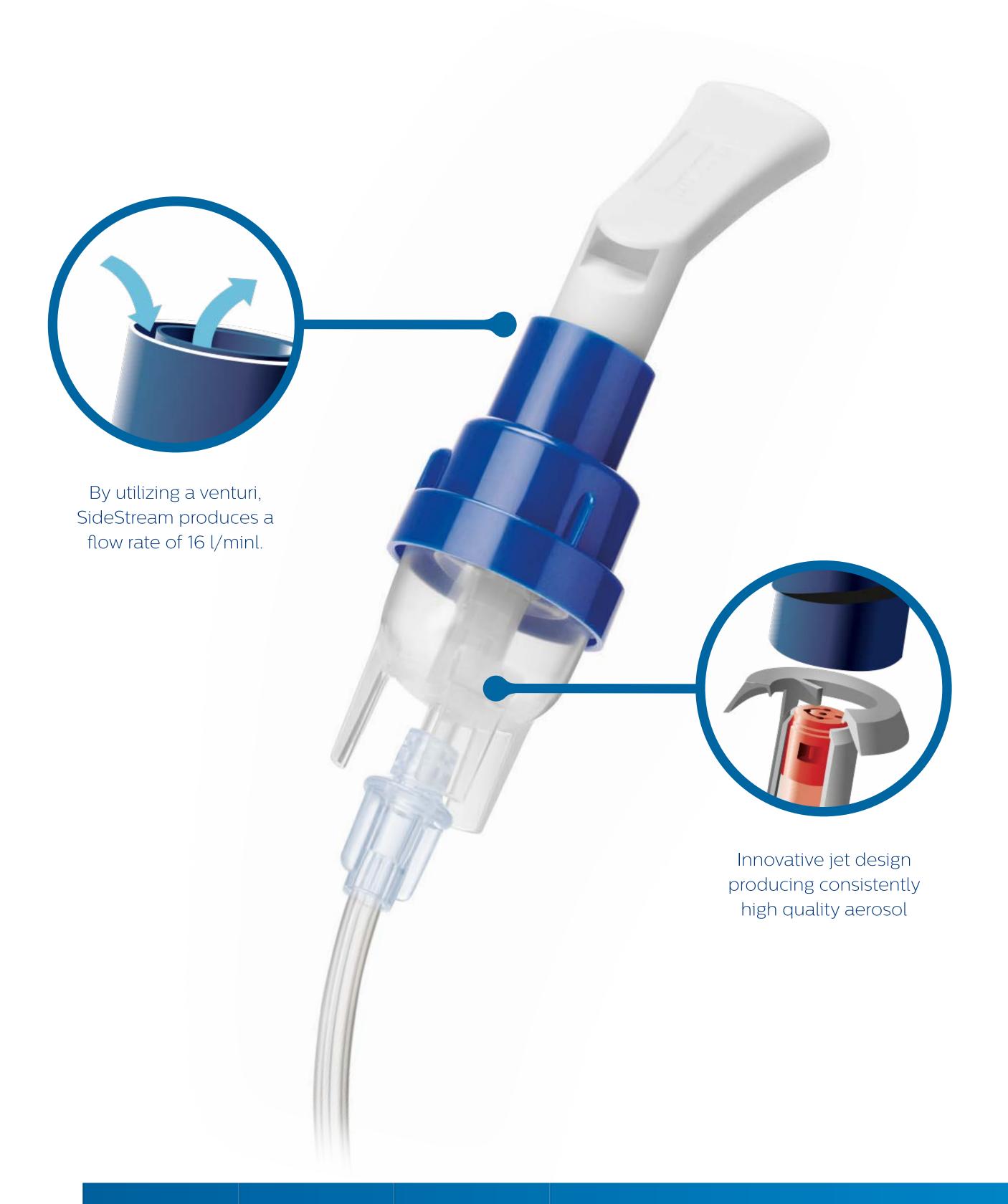

















Product specification

Product

Dimensions	6.7x9.3x3.8 inches
Weight	2.6lbs (1.17 Kg)
Noise level	60 +/- 4dBA
Mode of operation	Intermittent (30 min ON/30 min OFF)

External packaging weights and dimensions

Product description	Size	Weight
Philips Home Nebulizer (unit carton) 1 EA	8.8x5.8x6.5 inches	3.32lbs (1.51 Kg)
Philips Home Nebulizer shipped in cases of 6	18.5x9.4x14.2 inches	21.1lbs (9.6 Kg)















Ordering information

120V/60Hz US plug

Product description	Reorder number	UOM	NDC number	HCPCS code
Philips Home Nebulizer with SideStream reusable and disposable nebulizers and mouthpiece	1130530	1EA	08373-0530-00	E0570
Philips Home Nebulizer with SideStream disposable nebulizer and mouthpiece	1130531	1EA	08373-0531-00	E0570
SideStream reusable kit	HS860	1 case of 10	N/A	A7005
SideStream disposable kit	HS800	1 case of 50	N/A	A7003
Replacement filters	1083128	1 pack of 5	08373-0528-00	A7014
Optional accessories:				
Adult masks	1025531	1 case of 50	N/A	A7015
Pediatric masks	1025529	1 case of 50	N/A	A7015













What's included...

- 1 Home Nebulizer Compressor
- 2 SideStream Nebulizer
- 3 Mouthpiece
- 4 Nebulizer tubing
- 5 Filter

Optional accessories

- Adult and pediatric aerosol masks
- Replacement filters (sold as packs of 5)















References

- 2.5ml Salbutamol solution drug
 An evaluation of two aerosol delivery systems for rhDNase PL Shah et al. European Respiratory Journal 1997

